


ulti med Influenza A/B Test Cassette - Comparison of Flocked Swabs to Polyester Standard Swabs

The new swab in the ulti med Influenza A/B Card offers several advantages:

- Perpendicular Nylon fibers act like a soft brush and allow improved collection of cell samples.
- Capillary action between the Nylon fiber strands facilitates uptake of liquid samples.
- Sample stays close to the surface allowing easier elution.

Thus, the flocked swab increases the sensitivity of the Influenza A/B antigen test kit.

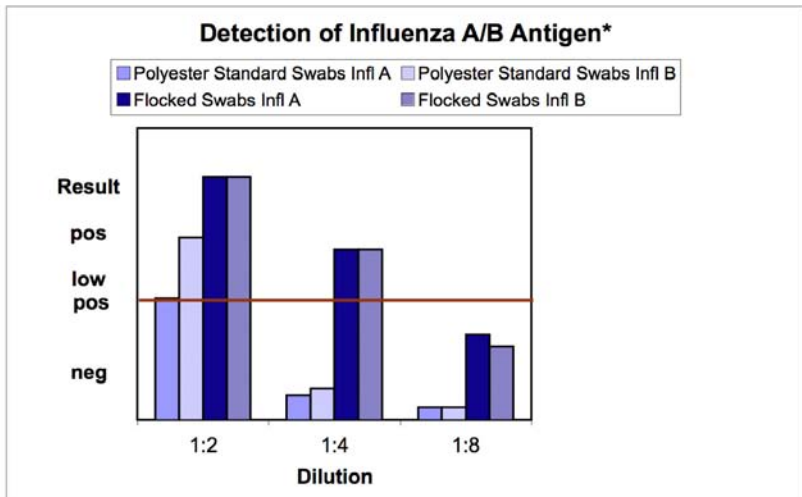

* The comparison of the swabs was performed with Influenza A/B test 004A474 and different dilutions of Influenza A and Influenza B positive strains