


## ulti med Streptococcus B Test Cassette - Comparison of Flocked Swabs to Polyester Standard Swabs

The new swab in the ulti med Streptococcus B Card offers several advantages:

- Perpendicular Nylon fibers act like a soft brush and allow improved collection of cell samples.
- Capillary action between the Nylon fiber strands facilitates uptake of liquid samples.
- Sample stays close to the surface allowing easier elution.

Thus, the flocked swab increases the sensitivity of the Streptococcus B antigen test kit.

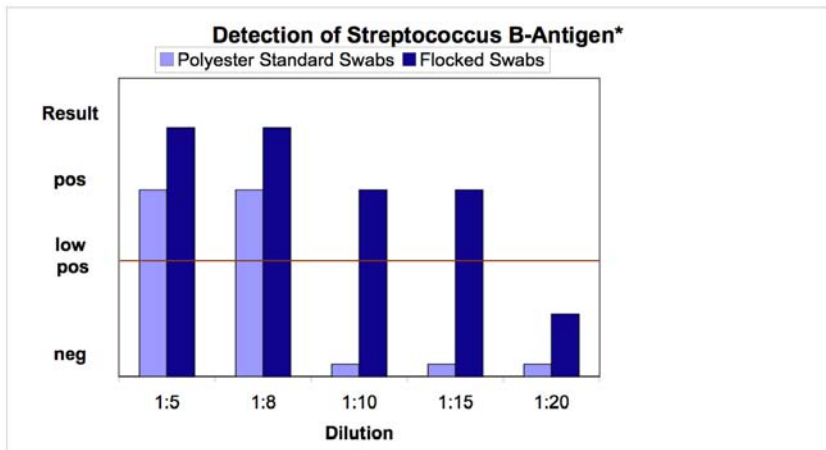

\* The comparison of the swabs was performed with ulti med Streptococcus B test # 014B280 and different dilutions of Strep B positive antigen control.