







## Benefits:

- 3.5" touch-screen display
- Optimized touch panel navigation
- Detection of up to 9 drugs in human saliva in one test
- Direct and fast evaluation
- Integrated backup function. Evaluation results will be stored in the device.
- GPS position tracking system 
- 2 USB connections 
- SD card reader (max. 32 GB) 
- Card reader for magnetic strips and micro chips (i.e. drivers licence)
- 4 x Sim Card Slots 
- RFID - Radio Frequency Identification Device / tag STD
- Fingerprint scanner (optional)
- Barcode reader (1D- and 2D-Barcode)

Able to read PDF 417 as used in drivers licenses

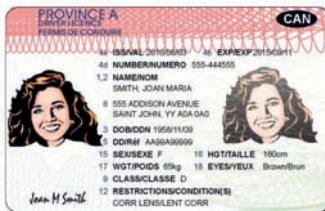

- GPRS internal (optional) 
- Leather protection bag
- Desktop application software

## Advantages

- reliable detection in saliva
- up to 9 drugs in just one test
- saliva samples almost impossible to manipulate
- quick and easy
- immediate positive results on THC consumption
- results will be stored on internal memory

## Our service for you

### Long time stability

In general, all ulmi med products are stable for a minimum of 18 months after production. Tests can be stored at room temperature, up to + 30°C.

### High reliability

The ulmi med rapid tests are extreme reliable and match in up to 99% (depending on test) with results of the reference method. They are also highly specific.

### Competent advice

Our team will be at your disposal for any kind of question concerning the handling of our products. Special questions concerning interaction and cross-reactivity with prescribed medicine cannot always be answered on the phone despite of our many years of experience. In this case please send your questions by fax or e-mail to our specialists. Usually you will receive a reply within 2 business days.

### Quality

ulmi med does not only fulfill the compulsory directions of their Quality Assurance System according to ISO EN 13485:2003, but also meets requirements which exceed these strict standards in particular with regards to extensive testing and complaint management.

### Further information

You will find all information concerning our products on our website [www.ulmimed.org](http://www.ulmimed.org) - of course, upon your request we will also send you any information by e-mail, fax or common mail.

Requests for product information please contact: [info@ulmimed.org](mailto:info@ulmimed.org)

### Research and Development / Production and Sales

Phone: +49 (0)4102-800 90  
Fax: +49 (0)4102-500 82  
e-mail: [info@ulmimed.org](mailto:info@ulmimed.org)  
Address: ulmi med Products  
(Deutschland) GmbH  
Reeshoop 1  
22926 Ahrensburg, Germany  
P.O. Box 18 14, 22908 Ahrensburg

### Sales

Phone: +32 (0)51 2004 25  
Fax: +32 (0)51 2004 49  
e-mail: [info@ulmimed.be](mailto:info@ulmimed.be)  
Address: ulmi med Products  
(Belgium) BVBA  
Honzebroekstraat 137  
8800 Roeselare, Belgium  
B.T.W. (V.A.T.) no. BE0451954672  
Register of Companies /  
HR Kortrijk 138687

### Quality assurance for your safety

ulmi med Products (Deutschland) GmbH complies with all requirements for production and sales of IVD (in vitro diagnostics) and is certified according to EN ISO 13485:2003.

Tax ID-Nr. DE 154 23 11 80

Register of Companies / Amtsgericht Ahrensburg: HRB 3562 founded 1992

Notified body TÜV SÜD Munich Germany

[www.drugreader.com](http://www.drugreader.com)

# SalivaScreen® Handheld Reader Touch

Digital saliva drug test analyzer,  
with integrated Printer and  
3.5" high-resolution touch-screen display



rev.24012012

Your first choice...

...a safe decision

Your first choice...

...a safe decision

## The SalivaScreen® Handheld Reader Touch

is a compact device that permits the digital interpretation and data management of saliva drug tests.

Therefore, it is the ideal device for roadside testing, emergency unit use, for emergency and corporate physicians, detention centres, border police etc.

**easy to use • secure • reliable**  
**complete with operating system**  
**optimized menu navigation with**  
**touch screen • integrated printer**



\* protection by utility patents

## The most important parts The Ultimate Test

This drug test device uses a unique cassette design, optimized for the most efficient detection of drugs in human saliva.

Saliva samples are easily accessible and can hardly be manipulated. Drugs can be detected after consumption even before they appear in a metabolized state in urine. Basically, the SalivaScreen Test has been designed with lower cut-offs as drug tests designed for the examination of urine samples. Saliva samples show a far better correlation with blood samples than urine samples in regards to their almost identical detection time frame.

### Up to 9 drugs out of 16 can be detected per test

Amphetamine (AMP), Benzodiazepine (BZO), Buprenorphine (BUP),  
Cocaine (COC), Cotinine (COT), EDDP (EDDP), Ketamine (KET),  
Marijuana (THC) 12 Metabolite, Marijuana (THC) 50 Parent,  
Methadone (MTD), Methamphetamine (MET), Opiates (OPI),  
Oxycodone (OXY), Phencyclidine (PCP), Barbiturate (BAR), Propoxyphene (PPX)

## The Software

The SalivaScreen evaluation software has been written and designed by german specialists and allows for

**intuitive handling • efficient workflow**  
**customized set-ups • fast results**

The software can be customized in it's own settings, but can also be tailor fitted to customer or authority requirements. The following gives a quick overview to the main drugtesting procedure.



Animations during the process are guiding the operator, but can be disabled for experienced users.



The Terminal Application is created to download the saved Testresult from the Handheld Device to the stationary computer.

The results are stored and displayed in a userfriendly way. The files are non editable and are for documentary purposes only. Future developments may include the connection to other databases.

Your first choice...

...a safe decision

Your first choice...

...a safe decision